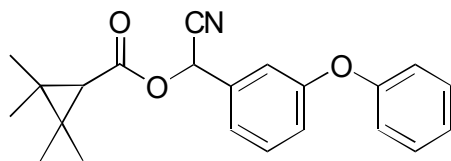


F-03

Fenpropathrin

**OA-4600** $\alpha=1.22$

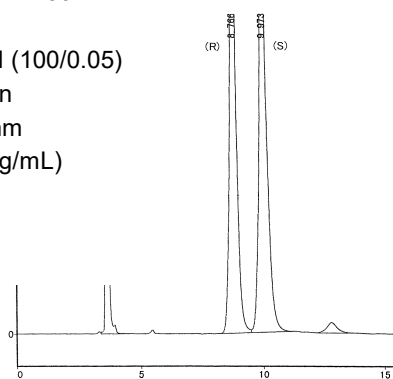
Column : 4.6mm i.d.×250mm

Mobile phase:

hexane/2-propanol (100/0.05)

Flow rate : 1mL/min

Detector : UV 230nm

Injection : 1 μ L (1mg/mL)**OA-4000** $\alpha=1.12$

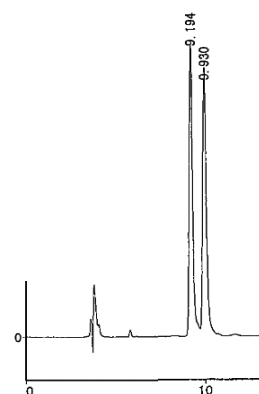
Column : 4.6mm i.d.×250mm

Mobile phase:

hexane/2-propanol (100/0.05)

Flow rate : 1mL/min

Detector : UV 254nm

**OA-4700** $\alpha=1.10$

Column : 4.6mm i.d.×250mm

Mobile phase:

hexane/2-propanol (100/0.05)

Flow rate : 1mL/min

Detector : UV 230nm

Injection : 1 μ L (1mg/mL)