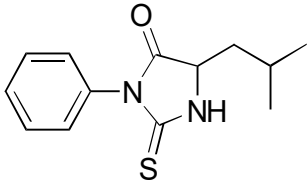
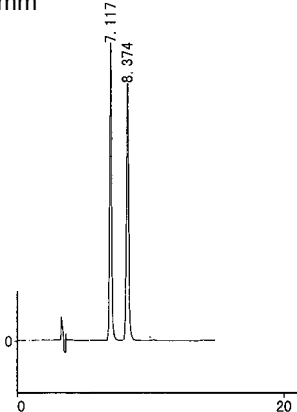


C-19	PTH-leucine
	
<p>OA-4700 $\alpha=1.33$</p> <p>Column : 4.6mm i.d.×250mm Mobile phase : hexane/2-propanol/ methanol/trifluoroacetic acid (90/8/2/0.2) Flow rate : 1mL/min Detector : UV 254nm Injection : 1μL (1mg/mL)</p> 	<p>OA-4100 $\alpha=1.25$</p> <p>Column : 4.6mm i.d.×250mm Mobile phase : hexane/2-propanol/ methanol/trifluoroacetic acid (90/8/2/0.2) Flow rate : 1mL/min Detector : UV 254nm Injection : 1μL (1mg/mL)</p> 