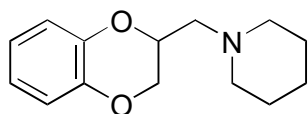


C-07

Piperoxan

**OA-4800** $\alpha = 1.31$ 


Column : 4.6mm i.d.×250mm

Mobile phase :

hexane/tetrahydrofuran/  
methanol/trifluoroacetic  
acid (60/30/10/0.2)

Flow rate : 1mL/min

Detector : UV 254nm

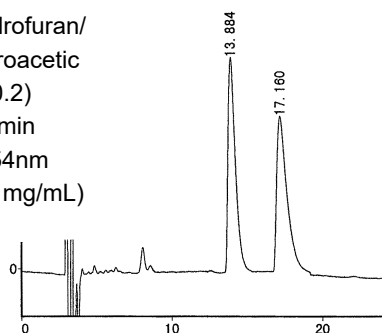Injection : 5 $\mu$ L (1mg/mL)**OA-4400** $\alpha = 1.23$ 

Column : 4.6mm i.d.×250mm

Mobile phase :

hexane/tetrahydrofuran/  
methanol/trifluoroacetic  
acid (60/35/5/0.2)

Flow rate : 1mL/min

Detector : UV 254nm

Injection : 5 $\mu$ L (1mg/mL)