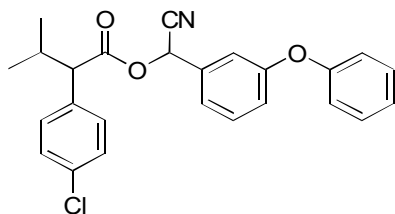


F-06

Fenvalerate

**OA-2000**

Column : 4.6mm i.d.×250mm


Mobile phase :

hexane/1,2-dichloroethane/
ethanol (500/30/0.15)

Flow rate : 1mL/min

Detector : UV 230nm

Injection : 1μL (5mg/mL)

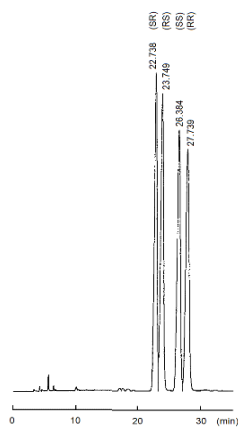**OA-2500I**

Column : 4.6mm i.d.×250mm

Mobile phase :

hexane/1,2-dichloroethane/
ethanol (500/20/0.1)

Flow rate : 1mL/min

Detector : UV 230nm

