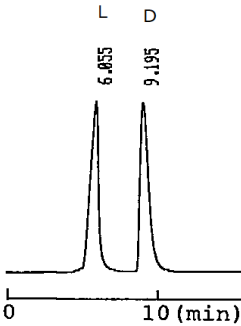


G-07	Aspartic acid	
<chem>NC(CC(=O)O)C(=O)O</chem>		
<p><b>OA-6000</b> <math>\alpha=1.68</math></p> <p>Column : 4.6mm i.d.×150mm Mobile phase : 1mmol/L copper( II ) sulfate in water Flow rate : 1mL/min Detector : UV 254nm</p> 	<p><b>OA-5000</b> <math>\alpha=1.37</math></p> <p>Column : 4.6mm i.d.×150mm Mobile phase : 2mmol/L copper( II ) sulfate in [water/2-propanol (95/5)] Flow rate : 1mL/min Detector : UV 254nm Injection : 5μL (5mg/mL)</p> 