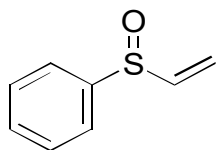


J-02

Phenyl vinyl sulfoxide

**OA-7600** $\alpha=1.18$


Column : 4.6mm i.d.×250mm

Mobile phase :

methanol/water (30/70)

Flow rate : 0.5mL/min

Detector : UV 254nm

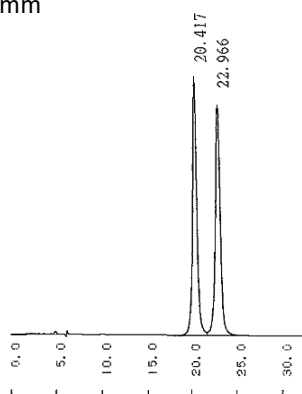Injection : 2 μ L (1mg/mL)**OA-2500I** $\alpha=1.03$

Column : 4.6mm i.d.×250mm

Mobile phase :

hexane/ethanol (95/5)

Flow rate : 1mL/min

Detector : UV 254nm

Injection : 1 μ L (1mg/mL)