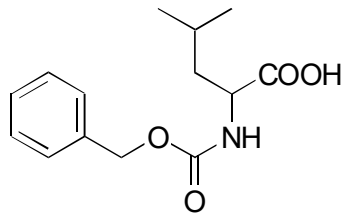


A-11

Z-leucine

**OA-3200** $\alpha=1.21$ 


Column : 4.6mm i.d.×250mm

Mobile phase :

0.01mol/L ammonium  
acetate in methanol

Flow rate : 1mL/min

Detector : UV 254nm

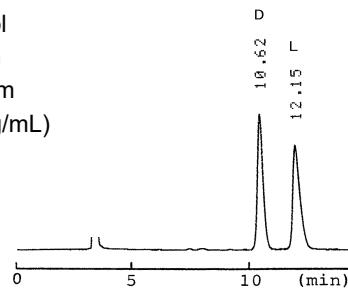Injection : 5 $\mu$ L (5mg/mL)**OA-3300** $\alpha=1.07$ 

Column : 4mm i.d.×250mm

Mobile phase :

0.01mol/L ammonium  
acetate in methanol

Flow rate : 1mL/min

Detector : UV 254nm

