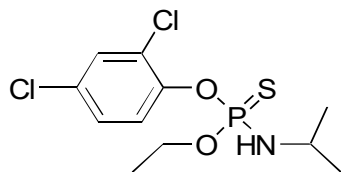


J-01

O-Ethyl O-2,4-dichlorophenyl isopropylphosphoamidethioate

**OA-4700** $\alpha=1.23$


Column : 4.6mm i.d.×250mm

Mobile phase:

hexane/2-propanol (100/0.05)

Flow rate : 1mL/min

Detector : UV 254nm

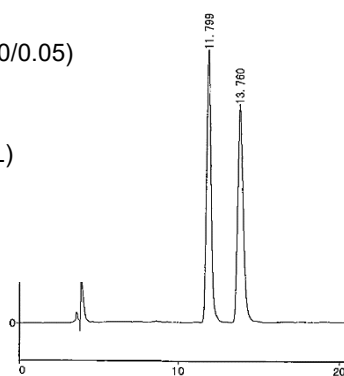Injection : 2 μ L (1mg/mL)**OA-4100** $\alpha=1.17$

Column : 4.6mm i.d.×250mm

Mobile phase:

hexane/2-propanol (100/0.05)

Flow rate : 1mL/min

Detector : UV 254nm

Injection : 2 μ L (1mg/mL)